



Single-use Ureteral Access Sheath

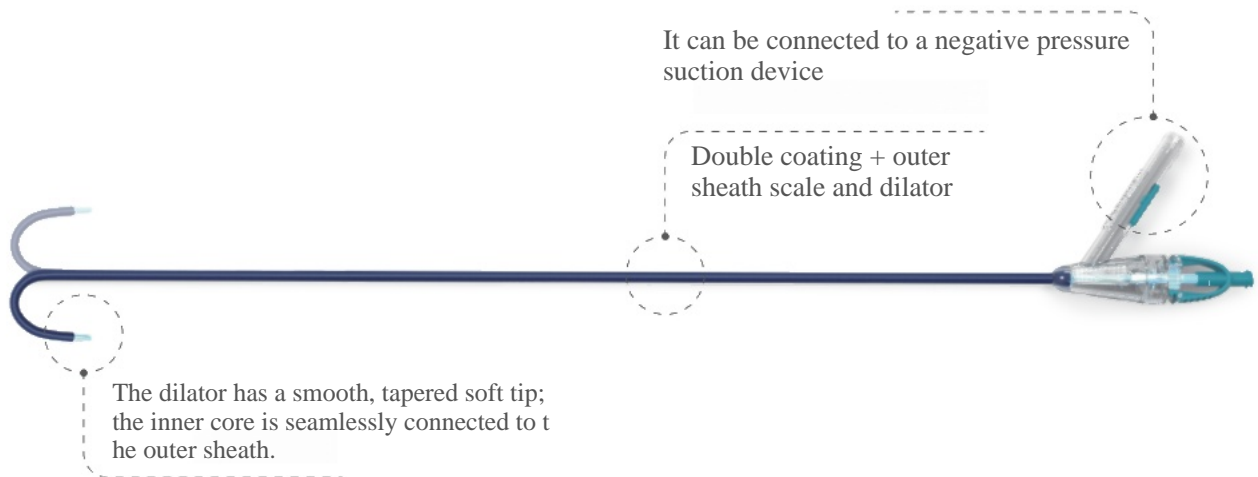

Demonstration with flexible ureteroscope

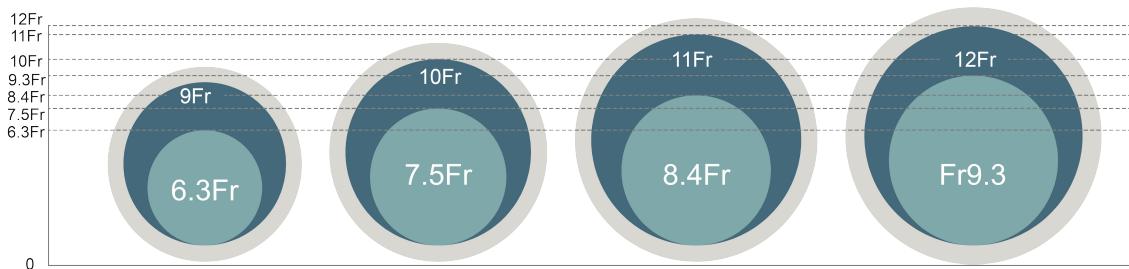

Product Feature


1. Both dilator and the guide sheath are coated with a hydrophilic coating, easy to insert and expand.
2. It can be connected to a negative pressure suction device, which can effectively reduce the intrarenal pressure, make the surgical field clearer, shorten the operation time, and improve the stone removal rate.
3. The spring-reinforced structure of the catheter body has anti-torsion and anti-compression properties, which ensures that the catheter passes smoothly and powerfully through the tortuous anatomical structure, and the inner wall is coated with a PTFE lubricant layer to protect the endoscope from damage ect.




Product Parameter

	Working length (mm)	Inner diameter (mm)	Obturator length (mm)
GYUAS-1036	360	3.4	450
GYUAS-1046	460	3.4	550
GYUAS-1136	360	3.85	450
GYUAS-1146	460	3.85	550
GYUAS-1236	360	4.07	450
GYUAS-1246	460	4.07	550
GYUAS-1336	360	4.4	450
GYUAS-1346	460	4.4	550
GYUAS-1436	360	4.7	450
GYUAS-1446	460	4.7	550

Vaina de Acceso Ureteral de un Solo Uso

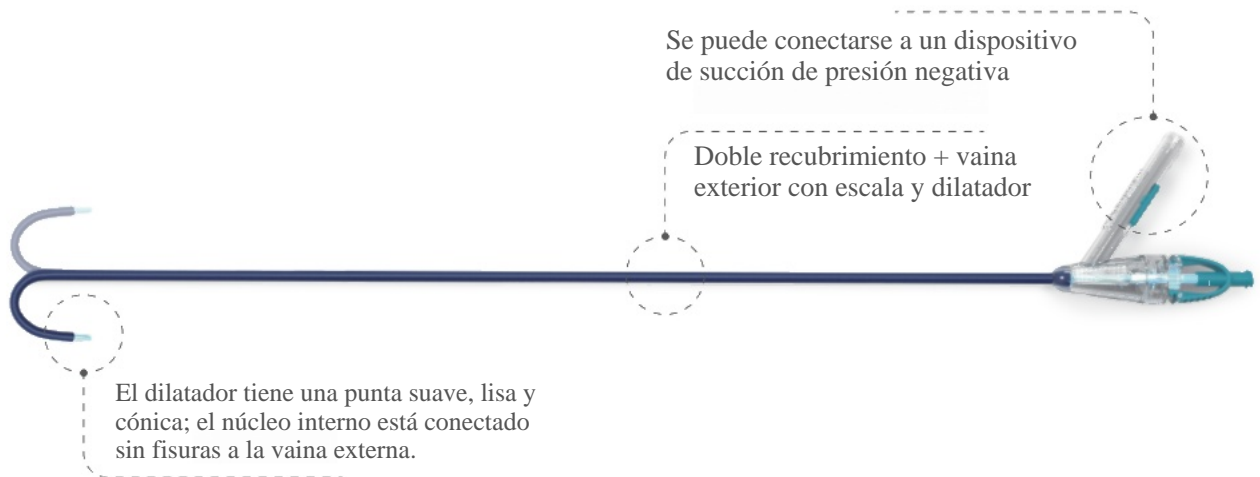

Demostración con ureteroscopio flexible

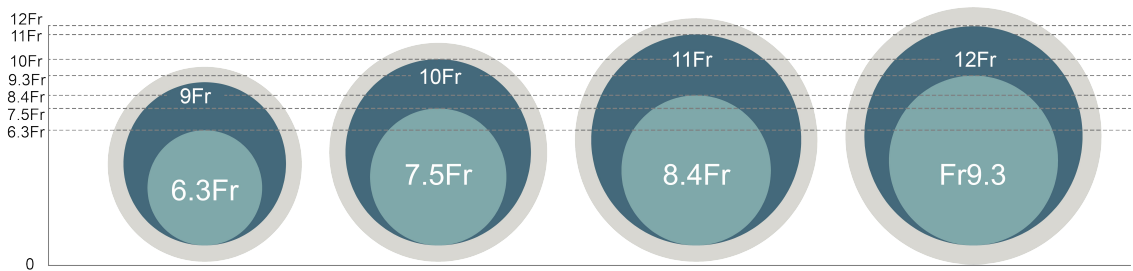

Característica del Producto

1. Tanto el dilatador como la vaina guía están recubiertos con una capa hidrofóbica, lo que facilita su inserción y expansión.
2. Se puede conectar a un dispositivo de succión de presión negativa, que puede reducir eficazmente la presión intrarrenal, hacer que el campo quirúrgico sea más claro, acortar el tiempo de operación y mejorar la tasa de eliminación de cálculos.
3. La estructura reforzada con resortes del cuerpo del catéter posee propiedades antitorsión y anticompresión, lo que garantiza que el catéter pase de forma suave y potente a través de la tortuosa estructura anatómica, y la pared interior está recubierta con una capa lubricante de PTFE para proteger el endoscopio de daños, etc.



Parámetros del Producto

Modelos	Longitud de trabajo	Diámetro interior	Longitud del obturador
GYUAS-1036	360	3.4	450
GYUAS-1046	460	3.4	550
GYUAS-1136	360	3.85	450
GYUAS-1146	460	3.85	550
GYUAS-1236	360	4.07	450
GYUAS-1246	460	4.07	550
GYUAS-1336	360	4.4	450
GYUAS-1346	460	4.4	550
GYUAS-1436	360	4.7	450
GYUAS-1446	460	4.7	550